

BioCLIA Allergy Reagent Kit, G3

(Chemiluminescent Microparticle Immunoassay)

REF

MY01225 (25 Tests/kit)

MY01190 (50 Tests/kit)

MY01155 (100 Tests/kit)

INTENDED USE

BioCLIA Allergy Reagent Kit, G3 is an *in vitro* quantitative assay for the determination of specific IgE to Cocksfoot (*Dactylis glomerata*, G3) in human serum via chemiluminescent immunoassay method on the fully automatic BioCLIA instrument. It is intended for *in vitro* diagnostic use as an aid in the clinical diagnosis of IgE mediated allergic disorders specific to Cocksfoot (G3) in conjunction with other clinical findings, and is to be used in clinical laboratories.

SUMMARY AND EXPLANATION OF THE TEST

Cocksfoot grass pollen is a very common allergen inducing asthma, allergic rhinitis and allergic conjunctivitis. ^{1, 2, 3} Cocksfoot was reported to be an important aeroallergen in Spain, Greece, Turkey and Japan. ^{4, 5, 6, 7}

In the *Dactylis glomerata* plant, allergenic components are most highly concentrated in the pollen, but are also present in the leaves and show traces in the stems. ⁸ This grass contains Group 1 allergens, to which more than 95% of patients allergic to grass pollen possess IgE antibodies. These are highly cross-reactive glycoproteins exclusively expressed in the pollen of many grasses. ^{9, 10, 11}

Cocksfoot grass also contains a Group 4 allergen. Almost 75% of patients allergic to grass pollen display IgE reactivity to Group 4 allergens, which hence can be regarded as major grass pollen allergens. ¹² Inhibition studies of IgE antibody binding to Dac g 4 (*Dactylis glomerata* - Cocksfoot grass) and to other pollen extracts confirmed the presence of cross-reactive allergens in *Secale cereale* (Cultivated Rye), *Lolium perenne* (Rye grass), *Festuca elatior* (Meadow Fescue), *Holcus lanatus* (Velvet grass), *Bromus arvensis* (Field Brome), *Poa pratense* (Meadow grass), *Hordeum sativum* (Barley), and *Phleum pratense* (Timothy grass). ¹³

Cocksfoot grass pollen also contains a Group 5 allergen. Almost 90% of grass pollen-allergic patients are sensitized to Group 5 grass pollen allergens. Group 5 allergens have been detected in Timothy, Rye grass, Meadow grass and Cocksfoot extracts.

PRINCIPLE OF THE PROCEDURE

BioCLIA Allergy Reagent Kit is a two-step immunoassay. The allergen of interest is biotinylated, it reacts with the specific IgE in patient sample and form a IgE-allergen complex with streptavidin beads. After incubation, a washing step removes the unbound reagents and sample matrix. Subsequently add enzyme conjugated anti-human IgE antibodies and bind to the IgE-allergen complex. After incubation, a second washing step removes the excess enzyme conjugates. Then addition of substrate results in the emission of light and RLU is measured. RLU is proportional to the amount of allergen specific IgE present in the sample. To evaluate the test results, RLU is transformed to concentration with the use of a calibration curve.

WARNINGS AND PRECAUTIONS

- For professional *in vitro* diagnostic use only.
 - Used on BioCLIA 6500 and BioCLIA 500 instruments only.
 - Do not use reagents beyond the expiration dates.
 - The kit contains human sourced materials. All recommended precautions for the handling of blood derivatives should be taken. Please refer to the OSHA Standard on Blood borne Pathogens¹⁰ or Biosafety Level ^{2,11} or other appropriate biosafety practices. ^{12,13}
 - Liquid waste and solid waste taken from BioCLIA 6500 and BioCLIA 500 should also be handled in accordance with the appropriate biosafety practices same as stated above.
 - Once opened, the reagent cartridge must be stored in the instrument's reagent carousel. For the first placement of reagent into the instrument, please take care to avoid spilling the reagents.
 - Spilled reagents should be cleaned up immediately. Comply with all federal, state and local environmental regulations when disposing of wastes.
- Improper cleaning or rinsing of the instrument may lead to chemical contamination of the reagents. Residues from common laboratory chemicals such as formalin, bleach, ethanol, or detergent can cause interference in the assay. Be sure to follow the recommended cleaning procedure as outlined in the BioCLIA 6500 and BioCLIA 500 User's Manual.



Precautions

- The allergen contains ProClin300 0.1% as preservative and Triton X-100 0.01% as surfactant and may cause an allergic skin reaction by skin contact. Avoid contact with skin. Wear protective gloves, protective clothing and take eye protection.
- Any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the competent authority of the Member State in which the user and/or the patient is established.

MATERIALS SUPPLIED

- G3 Allergen** Biotinylated G3 antigen in 0.01 M PBS (pH7.4) buffer with stabilizer.

G3	BIOT
----	------

Preservatives: 0.1% ProClin300;

Surfactant: 0.01% Triton X-100.

Kit size	Vial Quantity	Volume/Vial (mL)
25 Tests/Kit	1	0.375
50 Tests/Kit	1	0.75
100 Tests/Kit	2	0.75

ADDITIONAL MATERIALS SUPPLIED SEPARATELY

Product	CATALOGUE No.
BioCLIA Allergy Specific IgE Conjugate & Microparticle Set	MY01311 (100 Tests/Kit) MY01312 (150 Tests/Kit) MY01313 (200 Tests/Kit) MY01314 (300 Tests/Kit)
BioCLIA Allergy Specific IgE Calibrator and Control Set	MY01315
BioCLIA System Wash Buffer	MY00404
BioCLIA System Substrate	MY00405

BioCLIA 6500	MA00243
BioCLIA 500	MA00502
BioCLIA Cuvettes	MA00244 (2000 pcs/bag) MA00549 (65 pcs/box)
BioCLIA Silicone gasket (Small)	MV00195
BioCLIA Silicone gasket (Large)	MV00196

MATERIALS REQUIRED

- Distilled or deionized Water

STORAGE AND STABILITY

- Store the kit at 2-8 °C.
- The shelf life of the unopened kit is 18 months from day of production.
- Vial opened reagents or onboard reagents could be used for 28 successive days. The software of the BioCLIA instruments monitors the onboard (in-use) expiration of the reagent cartridge. The system will not allow use of a reagent which has expired.

SPECIMEN COLLECTION, STORAGE AND HANDLING

- Serum from venous can be used.
- Collect blood specimens using standard procedures.
- Test serum should be clear and non-hemolysis.
- Cloudy samples should be clarified by 3000r/min for 10 minutes before use. For samples with the presence of fibrin, ensure that complete clot formation has taken place prior to centrifugation of samples. Some samples, particularly those from patients receiving anticoagulant therapy, may require increased clotting time.
- Specimens may be refrigerated at 2-8 °C for up to seven days or stored at -20°C up to 36 months.
- Specimens may be kept onboard on BioCLIA instruments under room temperature (18-25°C) for up to 8 hours.
- Avoid repeated freezing and thawing.

ASSAY PROCEDURE

Detailed information about operating the BioCLIA instruments can be taken from the Instrument User's Manual.

Note that, it is important to perform all routine maintenance procedures for optimal performance.

Assay Calibration

The BioCLIA Allergy Reagent Kit utilizes seven points of calibrators to get the Working Calibration Curve. The calibrators information is uploaded into the instrument via the barcode provided in the BioCLIA Allergy Specific IgE Calibrator and Control Set. The Working Calibration Curve is generated and is used to calculate the analyte concentration (IU/mL) from the RLU obtained for each patient.

For each new lot of reagent, please calibrate prior to the first time use, and every 28 days thereafter. The software will not allow the lot to be used if the above requirements are not meet.

Programming and Running Samples

1. Put the reagents into the corresponding position of the reagent chamber of the fully automatic chemiluminescence analyzer. The information of the reagents can be uploaded into the instrument system through the scanning of reagent barcode, and can also be set through the supporting software of the

instrument.

2. The information of calibrator / quality control is identified by scanning the calibrator / control barcodes, and the position of calibrator / quality control is assigned in the instrument system.
3. The sample to be tested is placed on the instrument sample rack chamber, and the corresponding test information is edited through the instrument supporting software.
4. Start the operation procedure, and all calibrator / quality control / sample processing steps will be automatically executed.

Control

The control procedure could be done before running the specimens every day. Users also can adjust the control procedure period according to their own lab frequency.

CALCULATION OF RESULTS

Calculation and interpretation of results will be performed automatically by software on BioCLIA instruments.

RESULT INTERPRETATION

Classification System (Scoring):

Class	Conc. (IU/mL)	Interpretation
0	0-0.34	Negative/Equivocal
1	0.35-0.69	Low
2	0.70-3.49	Moderate
3	3.50-17.49	High
4	17.50-49.99	Very High
5	50.00-99.99	Very High
6	≥ 100	Very High

The test results only reflect the amount of allergen specific IgE presented in the collected sample and should be diagnosed in conjunction with other laboratory and clinical findings.

CUT-OFF VALUE DETERMINATION

Clinical significance cut-off of allergen-specific IgE concentration has been set at 0.35 IU/mL for indication for sensitization. BioCLIA allergy reagent kits have been set using this as the cut-off to differentiate positive versus negative results.

PERFORMANCE CHARACTERISTICS

APPEARANCE

Kit components are complete with no liquid leakage, no precipitation or floc in liquid reagents. Packing labels should be clear and easy to identify.

ACCURACY / SPIKED RECOVERY

The assay is designed to have a mean recovery within 100±15% when analyzing samples spiked with known amounts of allergen specific IgE. A positive serum with high value was spiked into a sample matrix (with low allergen specific IgE concentration, around 0.35 IU/mL) at the volume ratio of 1:9 and been tested in triplicates. Recovery rate was obtained by calculating the ratio between measured concentration to theoretical concentration. *

Sample	Obs-1	Obs-2	Obs-3	Av. Obs	CV	Expected	Obs/Exp
High positive serum	64.75	63.99	66.60	65.11	2.1%		

Sample matrix	0.92	0.89	0.84	0.88	4.8%		
Spiked sample	7.46	7.37	7.02	7.29	3.2%	7.31	99.7%

*Representative data; results in individual laboratories may vary from these data.

TRACEABILITY

The IgE calibrators can be traced to WHO International Standard Immunoglobulin E (IgE), human serum (NIBSC code: 11/234).

PRECISION

A study based on guidance from (NCCLS) document EP-A¹⁸ was performed.

Intra-assay precision: Four samples (negative, low, moderate, and high) were taken and tested with 10 replicates for each in a single run. Coefficient of variation (CV) was calculated for each of four samples. The results for intra-assay precision are shown in the table below.

Inter-assay precision: Four samples (negative, low, moderate, and high) were taken and tested with 4 replicates in a single run, two runs per day for 10 days. Coefficient of variation (CV) was calculated for each of four samples. The results for inter-assay precision are shown in the table below. *

Intra-assay precision: CV ≤ 10%				
Intra-Assay	negative	low	mid	high
Mean	0.31	3.59	18.91	81.62
CV	5.8%	2.7%	7.0%	2.0%

Inter-assay precision: CV ≤ 15%				
Inter-Assay	negative	low	mid	high
Mean	0.30	3.66	19.89	83.74
CV	9.7%	7.9%	11.1%	4.0%

*Representative data; results in individual laboratories may vary from these data.

LIMIT OF BLANK / DETECTION (LoB/LoD)

LOB/LOD was determined according to CLSI EP17-A guideline.

The assay is designed to have LoB/LoD of ≤ 0.1 IU/mL.

LIMIT OF QUANTITATION (LoQ)

The assay is designed to have LoQ at 0.1 IU/mL.

Linearity

The linear range of the assay is 0.1 - 100 IU/mL.

The linear range was determined by serially diluting a sample containing high levels of allergen specific IgE with a negative sample and covering the entire assay linear range according to the scheme in CLSI EP6-A. The expected value was plotted against the observed value, and linear regression analysis was performed to get slope, intercept and coefficient of correlation (r) values. The results are summarized in the table below*:

Slope	Intercept	r
1.01	-0.63	0.99

*Representative data; results in individual laboratories may vary from these data.

INTERFERENCE

No interference has been observed with bilirubin, hemoglobin, triglycerides, rheumatoid factor (RF), human anti-mouse antibody

(HAMA) at the level indicated below: *

- Bilirubin ≤ 200 mg/dL;
- Hemoglobin ≤ 400 mg/dL;
- Triglycerides ≤ 2,000 mg/dL;
- Rheumatoid factor (RF) ≤ 1,000 IU/mL;
- Human anti-mouse antibody (HAMA) ≤ 2,000 ng/mL.

CROSS REACTIVITY

One allergen specific IgE blank sample was spiked with the IgA/IgG/IgM to a different concentration and tested with the assay. No false positive results are observed for the samples spiked with IgA/IgG/IgM at 700,7000,500µg/mL, respectively.

INHIBITION STUDY

A study based on guidance from CLSI I/LA 20-A2 was performed.

Other allergens (C6, D2, E6, F92, I4, K202, M2, T6, W21) were added to one G3 positive serum at high concentrations (≥100 times of the G3 Allergen working concentration) to see if they interfere the performance of the assay. A recovery percentage for other allergens extract mixed samples should within 1±15%. The results for inhibition study are shown in the table below. *

Allergens in other categories	Spiked extract Conc. (ug/mL)	Inhibition rate
C6 extract	100	-12.1%
D2 extract	100	-9.9%
E6 extract	100	-1.9%
F92 extract	100	-2.3%
I4 extract	100	2.3%
K202 extract	100	-2.3%
M2 extract	100	-3.8%
T6 extract	100	1.2%
W21 extract	100	0.4%

*Representative data; results in individual laboratories may vary from these data.

METHOD COMPARISON

Method comparison was implemented by comparing clinical sample results of the assay to the results of predicated assay. The results are shown in the table below. *

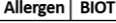
Clinical Sample	N=100	
Sensitivity	32/35	91.43%
Specificity	63/65	96.92%
Total agreement	94.18%	

*Representative data; results in individual laboratories may vary from these data.

LIMITATIONS

- The effectiveness of this kit is only confirmed for human sera, the applicability of the other kinds of samples is not verified.
- Reliable and reproducible results will be obtained when the assay procedure is carried out in accordance with the instructions and with adherence to good laboratory practice.
- Clinical diagnosis should not be made on the findings of a single test result, but should be interpreted with all clinical and laboratory findings.

SYMBOLS

	Catalog Number		Use-by date
	In Vitro diagnostic medical device		Lot Number
	Store between +2°C and +8°C		Consult Instruction for Use
	Manufacturer		Authorized Representative in the European Community
	CE Marking		Contains Sufficient for <n> Tests
	Biological Risk		GH507 Warning
	Biotinylated allergen		

REFERENCES

- Kitao Y, Sadanaga Y, Uno M, On N, Masuyama K, Ishikawa T. [Investigation of the genetics in allergic rhinitis. The 4th report--HLA class I and II specificities of orchard grass pollinosis]. Nihon Jibiinkoka Gakkai Kaiho 1988;91:516-20.
- Shimada T. [Four years study on Japanese cedar, orchard grass and ragweed pollinosis in Yotsukaido City--radioallergosorbent test (RAST) results of 361 patients]. Nippon Jibiinkoka Gakkai Kaiho 1986;89:864-71.
- Yamamoto Y, Tada R, Sakashita M, Sasabe T, Nakagawa Y, Yamada M, et al. [Allergic reactivity in cases of Japanese cedar and orchard grass pollenosis of the conjunctiva]. Nippon Ganka Gakkai Zasshi 1984;88:473-7.
- Subiza J, Cabrera M, Valdivieso R, Subiza JL, Jerez M, Jiménez JA, et al. Seasonal asthma caused by airborne Platanus pollen. Clinical & Experimental Allergy 1994;24:1123-9.
- Kontothanasi G, Moschovakis E, Tararas V, Delis A, Anagnostou E. Determination of sensitivity of inhalant allergens in patients with allergic rhinitis in West Athens. Rhinology 1995;33:234-5.
- Güneser S, Atici A, Cengizler I, Alparslan N. Inhalant allergens: as a cause of respiratory allergy in east Mediterranean area, Turkey. Allergologia Et Immunopathologia 1996;24:116-19.
- Shida T, Akiyama K, Hasegawa M, Maeda Y, Taniguchi M, Mori A, et al. [Change in skin reactivity to common allergens in allergic patients over a 30-year period. Association with aeroallergen load]. Arerugi = [Allergy] 2000;49:1074-86.
- D'Amato G, De PR, Verga A, Martucci P, Liccardi G, Lobefalo G. Antigenic activity of nonpollen parts (leaves and stems) of allergenic plants (Parietaria judaica and Dactylis glomerata). Annals of Allergy 1991;67:421-4.
- Grobe K, Becker WM, Schlaak M, Petersen A. Grass group I allergens (beta-expansins) are novel, papain-related proteinases. European Journal of Biochemistry 1999;263:33-40.
- Schenk S, Breiteneder H, Susani M, Najafian N, Laffer S, Duchêne M, et al. T Cell Epitopes of Phl p 1, Major Pollen Allergen of Timothy Grass (Phleum pratense). Advances in Experimental Medicine & Biology 1996;409:141-46.
- Hiller KM, Esch RE, Klapper DG. Mapping of an allergenically important determinant of grass group I allergens. Journal of Allergy & Clinical Immunology 1997;100:335-40.
- Fischer S, Grote M, Fahlbusch B, Muller WD, Kraft D, Valenta R. Characterization of Phl p 4, a major timothy grass (Phleum pratense) pollen allergen. Journal of Allergy & Clinical Immunology 1996;98:189-98.
- Leduc-Brodard V, Inacio F, Jaquinod M, Forest E, David B, Peltre G. Characterization of Dac g 4, a major basic allergen from Dactylis glomerata pollen. Journal of Allergy & Clinical Immunology 1996;98:1065-72.
- US Department of Labor, Occupational Safety and Health Administration, 29 CFR Part 1910.1030, Occupational Exposure to Bloodborne Pathogens.
- US Department of Health and Human Services. Biosafety in Microbiological and Biomedical Laboratories, Fourth Edition. Washington, DC: US

Government Printing Office, May 1999.

16. World Health Organization. Laboratory Biosafety Manual. Geneva: World Health Organization; 2004.

17. Clinical and Laboratory Standards Institute. Protection of Laboratory Workers from Occupationally Acquired Infections: Approved Guideline - Third Edition. CLSI Document M29-A3. Wayne, PA: Clinical and Laboratory Standards Institute, 2005.



HOB Biotech Group Corp., Ltd.
C6 Building, No. 218 Xinghu Road, Suzhou Industrial Park,
Suzhou, Jiangsu, 215123, China

CONTACT INFORMATION:

TEL (+86)512-69561996

Fax (+86)512-62956652

WEBSITE: www.hob-biotech.com

CUSTOMER SERVICE: TEL (+86)4008601202



EUROPE REPRESENTATIVE: Emergo Europe

ADDRESS/LOCATION:

Prinsessegracht 20, 2514 AP The Hague, The Netherlands



The eIFU is available on Website:

<http://en.hob-biotech.com/usercenter/login.aspx>

TECHNICAL ASSISTANCE

For technical assistance, contact your National Distributor.

Date of issue: 24th May 2021

Version: A/0 (EN)